


B-0015-18-BR RockHard 2018-EN



1112

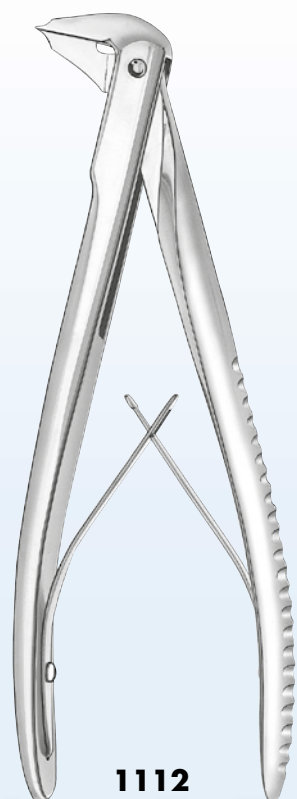

- smooth operating, precisely gripping crown remover
- once the crown has been cut through (size 0.10 - 0.12 crown cutter), the beak is inserted correctly (refer to sketch) and gentle pressure exerted to spread the crown. As the crown is spreaded parallelly, there is much less risk of fracturing the preparation.
- an Allen key comes with the crown remover for adjusting the beak if necessary (No. 1112 only).



only for
No. 1112



straight

LS783/18
18cm



curved

LS783A/18
18cm

LIQUIDSTEEL

Very fine MICRO tweezers with RH (Rock Hard) coating. Designed to grip the finest tissue atraumatically and to position membranes (without perforation)

ROCKHARD

the new coating series



LONG LIFE PERFORMANCE

by **ROCKHARD** the power coating

The diffusion coating made of Titanium Carbide with high hardness, stability and elasticity.


Experience the **hard-as-granite** needle holder!

The Titanium Carbide ions are deposited to the working end in a special process, which allows the ions being penetrated into the material.

As a result, the degree of hardness increases from normal 48-52 HRC (Rockwell) to 75-80 HRC. The wear of the surface is then reduced enormously. The risk of fractures which is possible on the conventional TC deposits, and bending which occurs especially on those very fine needle holders is eliminated almost completely.

All instruments from ROCKHARD series come with a blue handle.
(except No. 1112)

ROCKHARD Instruments

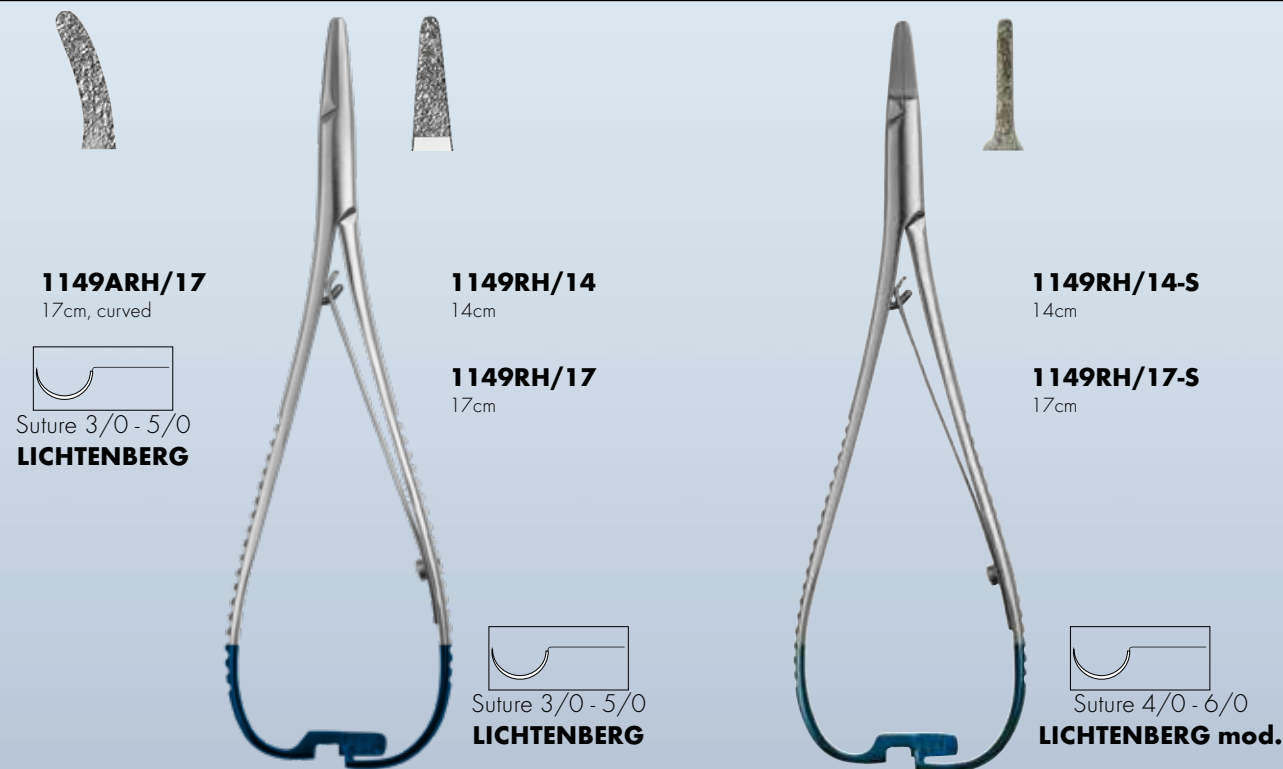

ROCKHARD Instruments

