





Die neuen Zangen mit BLISS!

Ein feiner, schlanker und seitlich abgeflachter Zangenkopf ermöglicht auch ein tiefes Fassen des Zahnes bis ins Knochenfach hinein, sowie eine optimierte, approximale Rotationsbewegung.

Die vom Zangenkopf distal geneigte Zahnung am Rand der Zangenmäuler erlaubt ein komfortables Vorschieben bis zum Zahnhals. Die neue Zahnungsstruktur bietet einen optimalen Halt für das Ausführen einer Rotationsbewegung sowie bei der Extraktion des Zahnes.

- medizinischer Edelstahl im CRYO-Härteverfahren –ice-tuned– veredelt
- elektrolytisch passiviert zur Vermeidung von Korrosion
- blendfreie Oberfläche durch spezielle Mattierung
- geeignet für Aufbereitung im Thermodesinfektor und Autoklav (134°C)

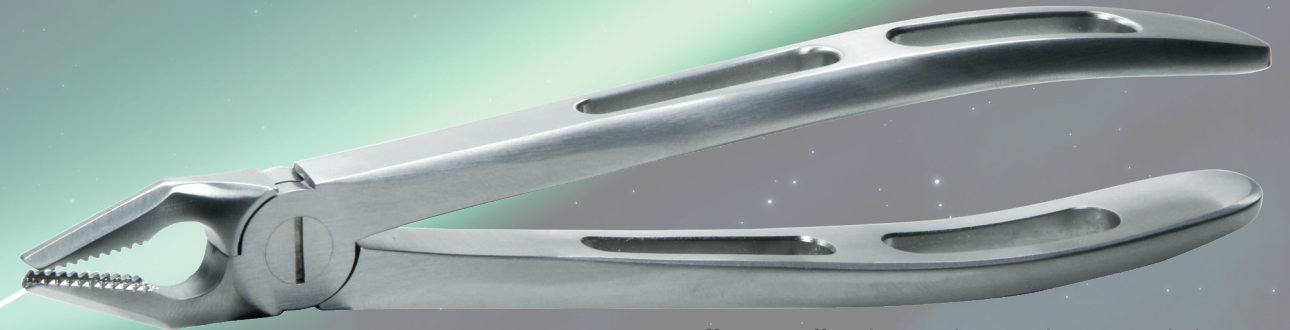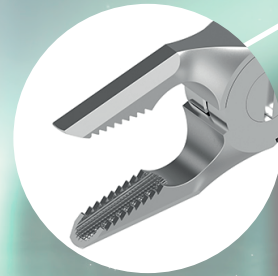

Die offene Griffstruktur reduziert das Gewicht* der Zange und ermöglicht eine gute Durchspülung insbesondere bei der thermischen Aufbereitung.
*(ca. 35% Gewichtsreduzierung gegenüber Standard-Zahnzangen aus unserem Programm)



GX1
für obere Schneide-
und Eckzähne



GX7
für obere Prämolaren




GX17
für obere Molaren
rechts



GX18
für obere Molaren
links




GX51A
für obere Wurzeln

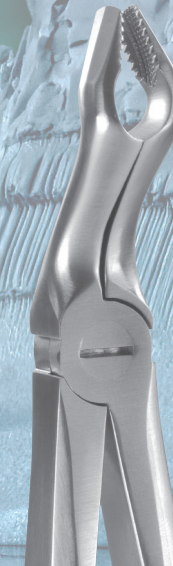

GX67A
für Weisheitszähne,
Oberkiefer



GX100
SET aus zehn
GatorEx®-Zangen
Oberkiefer-Fig. 1, 7, 17, 18, 51A, 67A
Unterkiefer-Fig. 13, 22, 33A, 79

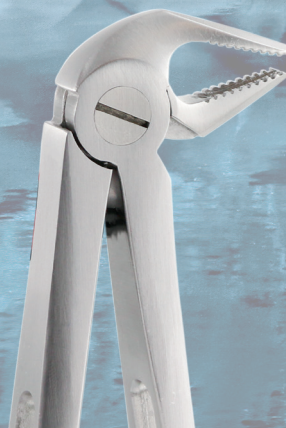

GX13
für untere Prämo-
la-
ren, Schneidezähne
und Eckzähne



GX22
für untere Molaren




GX33A
für untere Wurzeln



GX79
für Weisheitszähne,
Unterkiefer



GX136
für untere Wurzelreste

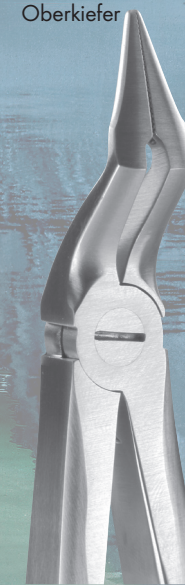

GX152
für obere Wurzelreste



GX86C
für untere Molaren
Kuhhornzange



FINEST DENTAL INSTRUMENTS
- SINCE 1916 -

NEUENKAMPER STR. 80 - 86 · D-42657 SOLINGEN · GERMANY
TEL. +49(0)212/81 00 44 · FAX +49(0)212/8 73 45
E-MAIL: INFO@CARLMARTIN.DE · INTERNET: WWW.CARLMARTIN.DE





Hier geht es zum Shop



The new pliers with BITE!

A fine, delicate and at the side flattened forceps head enables also deep grasping of the tooth in the bone pocket as well as an optimized, proximal rotation movement. The distally angled serration at the edge of the forceps' jaw enables a comfortable moving forward up to the neck of the tooth. The newly structured serration provides an optimal grip for executing rotational movements as well as during the extraction of the tooth. The open grip structure reduces the weight of the forceps and allows a thorough rinsing, particularly at the thermal processing.

- medical stainless steel refined in the CRYO hardening process – ice-tuned
- passivated electrolytically to avoid corrosion
- lare-free surface by special matting
- appropriated for processing in the thermal disinfectant and autoclave (134° C)

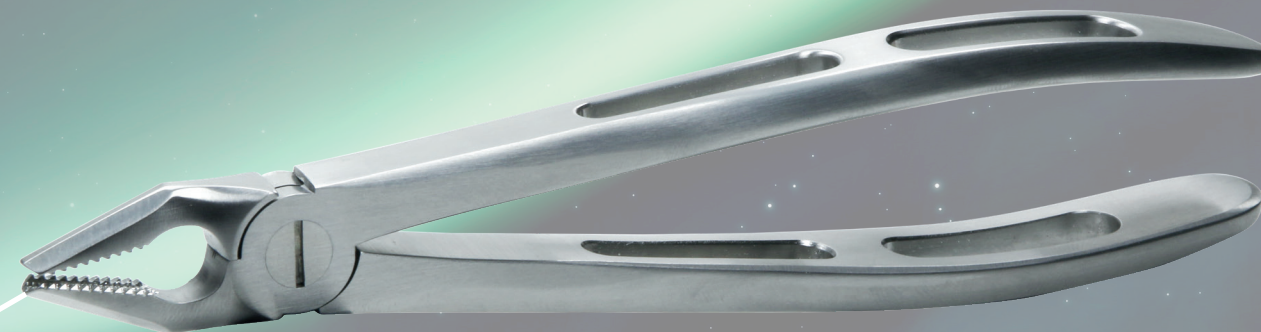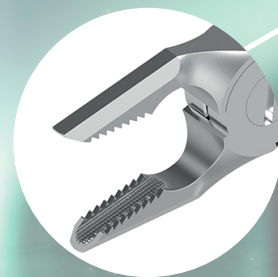


GX1
for upper incisors
and canines



GX7
for upper premolars




GX17
for upper molars
right

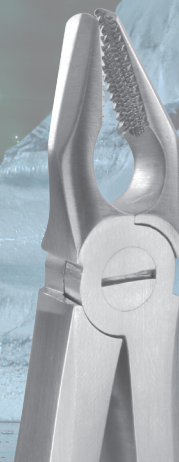

GX18
for upper molars
left



GX51A
for upper roots

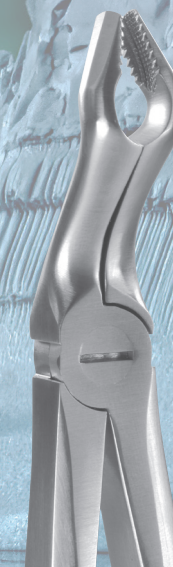

GX67A
for wisdom teeth,
upper jaw



GX100
SET of ten GatorEx® pliers
Upper jaw - Fig. 1, 7, 17, 18, 51A, 67A
Lower jaw - Fig. 13, 22, 33A, 79



GX13
for lower premolars,
incisors and canines



GX22
for lower molars



GX33A
for lower roots



GX79
for wisdom teeth,
lower jaw



GX136
for lower root fragments



GX152
for upper root fragments



GX86C
for lower molars
cow-horn beak



Shop now!



FINEST DENTAL INSTRUMENTS
- SINCE 1916 -

NEUENKAMPER STR. 80 - 86 · D-42657 SOLINGEN · GERMANY
TEL. +49(0)212/81 00 44 · FAX +49(0)212/8 73 45
E-MAIL: INFO@CARLMARTIN.DE · INTERNET: WWW.CARLMARTIN.DE